
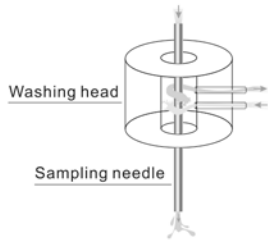
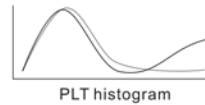
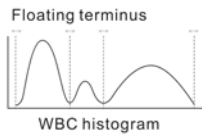


H-7022

3-PART 23 PARAMETERS HEMATOLOGY ANALYZER TOUCH SCREEN

H-7022 Hematology Analyzer is a 3-part differential of WBC blood cell counter. It is advantaged at digital technology, auto-relocation separatrix, intelligent flush clog, dual measuring technology (with time and volume) to ensure precise of the test results.


Advanced Technology

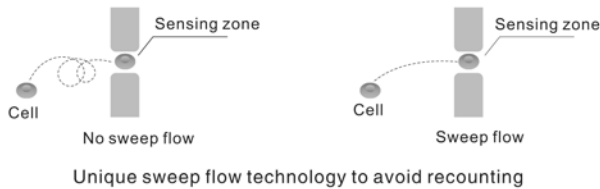

Auto sample probe cleaning



Auto standby and wakeup

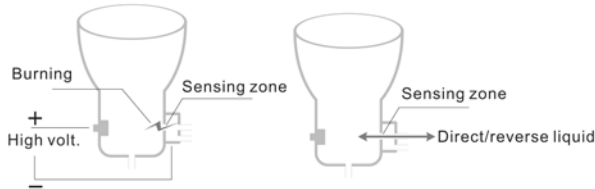

Automated & User-friendly Interface

- 8 inch LCD display, resolution:800*600
- Real-time alert for abnormal results
- USB keyboard & mouse operation
- Based on Linux
- Easy to search and edit patient information
- Storage capacity:50,000 sample results include histograms
- Optional printout for formats
- Real-time guideline for easy operation

Accurate Results & Low Sample Consumption

- Accurate analysis technology with low sample consumption:
 - Venous mode: 9.6ul
 - Capillary mode: 9.6ul
 - Pre-diluted mode: 20ul



Easy Maintenance & Environment-friendly System

- Adjustable sample probe for various cuvettes
- With liquid sensor to alarm low reagent
- Cyanide-free Lyse
- Automated & Systematic diagnosis and maintenance with alerts
- Liquid tubing and electronic circuit separately for easy maintenance
- Closed reagent tubing to avoid contamination

Perfect QC System

- Calibration mode: manual and auto
- QC mode: L-J, SD, CV, QC histogram

Three sample modes: venous, capillary and pre-diluted

23 parameters:



WBC

Lymph#

Mid#

Gran#

Lymph%

Mid%

Gran%

RBC

HGB

HCT

MCV

MCH

MCHC

RDW-CV

RDW-SD

PLT

MPV

PDW

PCT

P-LCR

